

# OroNasal 7450 V2 Mask

With Adapters From Hans Rudolph





OroNasal 7450 V2 Mask™ with Headgear and way T-shape™ Non Rebreathing Valve



OroNasal 7450 V2 Mask™ with Headgear and Custom Mask Adapter

## Product Features

- Five sizes of face pieces to cover a wide range of patient sizes.
- New facial sealing flange is a single piece mask made of clear soft silicone rubber.
- Anatomically contoured face mask design with resilient ribbed support, structure for a leak free comfortable fit and integrated chin cup for stability.
- Easily adjustable, light weight, 5 strap headgear for patient comfort, available in three sizes.
- New quick release strap clips for quick and simple mounting and dismounting of headgear.
- Reusable mask face pieces steam autoclaveable, cold chemical, pasteurisation and dishwasher safe.
- Reusable NRB valves validated for high level disinfection by cold chemical.
- NRB valves available in T & Y shape in several sizes covering a full patient range.
- Custom mask adapters designed and manufactured per your request

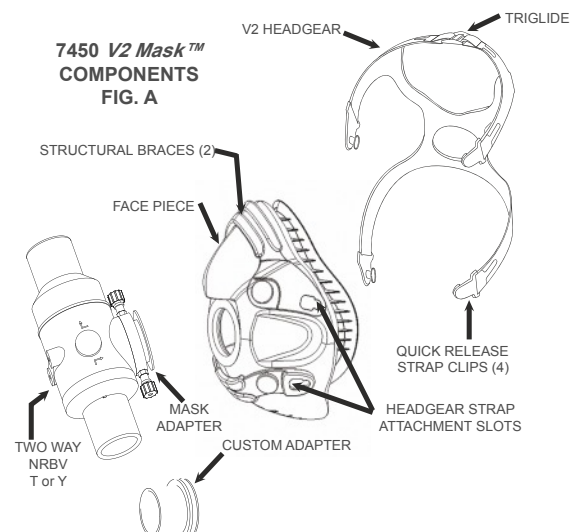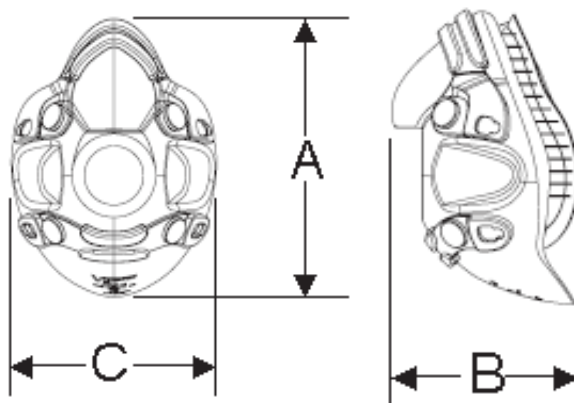

# Technical Specifications

## Mask and Headgear Service Life

The mask is expected to stay in service for minimum of 25 disinfection or sterilisation cycles or 6 months of use under normal conditions, whichever ever occurs first. The headgear is expected to stay in service for 6 months of use.

## Operational conditions for use

- Temperature Range: 540 C
- Humidity Range: 095% RH

## Cautions

Patients with facial hair, especially beards, may experience mask leakage which may require shaving hair. Do not use these valves in a life support system such as ventilator circuit.

## Materials

▪ Face Piece with Structural Braces
▪ Headgear
▪ Headgear Hook
▪ Headgear Strap Clips (4), Triglides
▪ Custom Mask Adapters
▪ Two Way Non Rebreathing Valves (NRBV)
▪ Silicone Rubber Blue, Polycarbonate Thermo Plastic
▪ Polyurethane Foam Black, Nylon UBL Gray & Nylon Fabric Red
▪ Nylon, Black
▪ Polypropylene, Black
▪ Polycarbonate Thermoplastic
▪ Refer to NRBV Data Sheet

## Face Mask Sizes

Characteristics	Large	Medium	Small	Extra Small	Petite
A (Inches)	5.9	5.6	5.2	4.8	4.6
B (Inches)	3.7	3.4	3.2	3.1	3.0
C (Inches)	4.2	4.2	4.0	3.9	3.9
Dead space (ml)	143	125	99	88	78
Weight (gm)	128	118	100	89	83

## Love Medical Ltd

Unit 1, Willan Enterprise Centre,  
Fourth Avenue, Village, Trafford Park,  
Manchester, M17 1DB, United Kingdom

www.lovemedical.com  
info@lovemedical.com  
sales@lovemedical.com

Tel: 0161 976 2744  
Fax: 0161 976 2764